

Experience the healing nature of light



MedX **DENTAL** Product
OraLase

Clinical Effects of Low Level Laser Therapy

- Muscle relaxation
- Desensitization
- Pain reduction after dental treatment
- Reduction of pain and swelling
- Internal wound healing
- Stimulation of fibroblasts, osteoblasts and odontoblasts
- Treatment of nerve injuries, included paraesthesia and trigeminal neuralgia
- Accelerate tissue repair
- Stimulation of wound healing

Benefits of MedX Phototherapy

- Treatment of large muscle groups
- Highly effective TMJ treatment
- Muscle relaxation at end of long procedures
- Continuous wavelength treatment
- Allows you to differentiate practice with effective and leading edge technology
- Effective in treating conditions that don't normally respond to traditional therapy interventions
- Simultaneous multiple site treatment options
- This healing method often eliminates the need for drug therapies which are usually accompanied by harmful side effects

Used For


- Endodontics
- Extractions and wound healing
- Facial pain
- Headache (tension or cervicogenic)
- Temporomandibular joint dysfunction (TMJ)
- Hypersensitivity
- Nerve damage
- Restorative procedures
- Implants
- Gingivitis and periodontitis
- Pre and post op pain



MedX low level laser therapy is a proven and highly effective treatment which reduces pain, decreases inflammation/edema, accelerates tissue repair and increases local circulation.



***** PH Toll Free 1-888-824-7558 *****
Visit Online <http://ColdLaserSupplies.com>

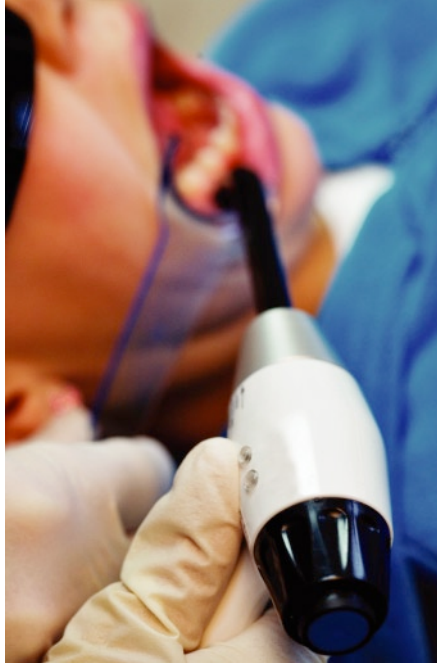

The MedX portable laser has one GaAlAs laser diode and a visible guide LED. The unit is lightweight with a built-in rechargeable Li-Ion battery.

Reduced pain and faster healing

MedX LPS200 Portable Dental Laser Specifications

CATEGORY	SPECIFICATION
Laser Diode	GaAlAs
Laser Equipment Safety Classification	Class IIIb
Nominal Laser Wavelength	808nm
Guide Light	Three 633nm visible guide LEDs
Duty Cycle	100% (Continuous mode)
Optical Output Power	200mW +/- 10% at 808nm
Max. Optical Output Power	250mW
Output	Continuous: 2J – 10 Seconds 4J – 20 Seconds 8J – 40 Seconds
Laser Power Monitoring	Automatic, continuous
Treatment Area at Skin	10 mm ²
Safety Interlock	Yes
Battery	3.6V DC Li-Ion, 1500 mAh capacity
Dimensions	23cm x 8cm x 4cm
Weight	140g (not including accessories)

The MedX Treatment Manual illustrates comprehensive protocols for a wide variety of common conditions, enabling the clinician to effectively treat patients, to accelerate healing and reduce pain for improved outcomes.

Phototherapy has gone through stringent regulatory processes. MedX products are US FDA and Health Canada cleared and are produced in an ISO 9001: 2000 and ISO 13485: 2003 certified manufacturing and testing facility.

MedX Portable Laser Kit Includes:

Hand-held low-level laser:

The MedX laser provides a total of 200mW of power with one GaAlAs laser diode and a visible guide LED. Class IIIb laser.

Intra-oral and Extra-oral Light Guides

Carrying case

MedX Treatment Manual:

Highly effective and efficient protocols illustrate how to deliver laser treatments for all dental conditions.

Operating manual

Two sets of goggles:

Patient and clinician goggles

Power supply for recharging unit:

9V DC Power Adapter



Distributed by

<http://ColdLaserSupplies.com>



***** PH Toll Free 1-888-824-7558 *****